







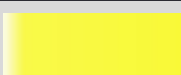





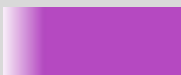
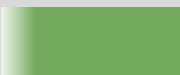







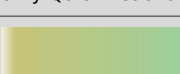

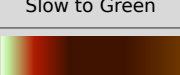
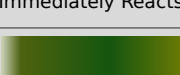
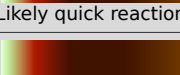
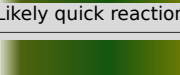
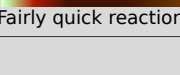


COLOR REACTIONS FOR MARQUIS, MANDELIN & MECKE REAGENTS

SUBSTANCE	MARQUIS	MANDELIN	MECKE	NOTES
MDMA, MDEA or MDA	 5sec or Less to React	 Immediately Reacts	 Immediately Reacts	>10sec to react with Marquis / Mandelin indicates very low potency
MDMA mixed with (Meth)amphetamine	 Immediately Reacts	 Immediately Reacts	 Immediately Reacts	May need to test a couple of times with Marquis to see yellow orange
(Meth)amphetamine	 Reacts Quickly	 Reacts Quickly	No Reaction	Methamphetamine and amphetamine reactions are the same with Marquis, Mandelin and Mecke
Methylone (bk-MDMA)	 Immediately Reacts	 Immediately Reacts	 Slowly Changes Colour	
Mephedrone (4-MMC)	No Reaction (Some impure batches have a yellow reaction)	 Immediately Reacts	No Reaction (may smoke)	
LSD	 No Reaction	Unknown (Likely no reaction)	Unknown (Possibly green) reaction	Use fairly large pill scraping or reaction won't be seen. Test won't work with LSD tabs (blotters).
Ketamine	No Reaction	 No Reaction	No Reaction	Sometimes PMA shows no reaction to Mecke and could be mistaken for Ketamine!
Opiates (such as Codeine / Morphine / Heroin)	 No Reaction	 No Reaction	 No Reaction	
2C-B	 Blue after a few mins	 Reacts quickly	 Reacts quickly	Check website for reactions to other 2C-x or 2C-T-x chemicals
2C-T-7	No Reaction	 Fairly Quick Reaction	 Fairly Quick Reaction	2C-T-7 -dangerous if consumed with MDxx* or a number of other drugs
DXM	 >10sec to react	 Slow to Green	 Immediately Reacts	DXM -dangerous if consumed with MDxx* or a number of other drugs
PMA	No Reaction	 Likely quick reaction	 Likely quick reaction	PMA - high risk of death if consumed
PMMA	No Reaction	 Fairly quick reaction	 Fairly quick reaction	PMMA - high risk of death if consumed
Piperazines (such as BZP / TFMPP / mCPP)	No Reaction	No Reaction	No Reaction	Piperazines - high risk of bad side effects if consumed and some risk of death
Strychnine	No Reaction	 No Reaction	No Reaction	Very toxic poison - high risk of death if consumed

EVEN USING THIS CHART THE REACTIONS MAY LOOK THE SAME FOR SOME DIFFERENT MIXTURES OF DRUGS AND OR IMPURITIES!

MDxx + Piperazines - will appear same reaction as pill containing only MDxx. High risk of bad side effects if consumed and some risk of death.

***MDxx - means MDMA or similar substance such as MDEA or MDA.**

CHART CREATED IN THE INTERESTS OF HARM REDUCTION. Please check pillreports.com and ecstasydata.org to see if there are any reports on your pill. Special thanks to Ekstasis for this chart.

Chart available for download at <http://testkitplus.com/drug-info/links/>